

# ORDERING PROCESS SERVICE+ “PRODUCTION”

Dental  
laboratory

Ivoclar Digital  
Service+




## 1. Sign in/registration

INBOX

3shape

Install the Ivoclar Digital CAD data base

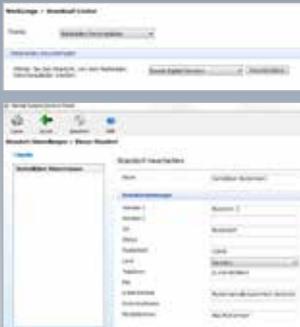

Additional CAD systems




Order Manager



## Ivoclar Digital material parameters

Software parameter charts are available for the following systems:

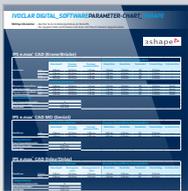

## Import files available for:

3shape

- Material parameters
- Implant systems
- Attachment/bar profile
- Connection to Inbox



- Material parameters
- Implant systems

exocad

- Implant systems
- Attachments
- For the material parameters please refer to the software parameter chart

More information at:

[www.ivoclarvivadent.com/en/cad-cam-download-center/](http://www.ivoclarvivadent.com/en/cad-cam-download-center/)



## 2. Scan & design

Individually design your restorations taking into account the respective material parameters



## 3. Data transfer/upload

### – Tooth-supported restorations



Select the manufacturer:  
"Ivoclar Digital Service+"  
and send the data



Create order and upload data

### – Implant-supported restorations (bars/bridges)

In order to produce implant-supported restorations, the master model is required.  
The implant positions will be measured in the Ivoclar Digital milling centre and compared with the scan data. A 100% fitting accuracy can therefore be guaranteed.

- Pack the model records securely (implant model incl. full, uncut gingival mask)



- Include the printout of the order form

### 3Shape



### Order Manager



### – Dispatch options

#### Pick-up order



#### Self-dispatch to

Ivoclar Vivadent Manufacturing GmbH  
Gustav-Flora-Straße 32  
39025 Naturns | BZ | Italy



## 4. Production

- Start of production with specified material
- Final quality check



## 5. Dispatching

UPS parcel dispatch to the laboratory